

VP Premium 2.0



Wearers demanding a personalized progressive lens that adapts to their way of looking and to their prescription.

Specially recommended for frequent users of hand held devices.



Lens personalized according to real position of wear



Adapted to the digital world



Optimized near vision

VP Premium 2.0

A visual solution for demanding wearers.



Personalized progressive lens
Frequent use of digital devices
Reading area improvement



4286 Mt. Carmel-Tobasco Rd., Suite E
Cincinnati, Ohio 45244 | **Tel:** 513-273-4034

www.RandOptical.com




R & D 
OPTICAL LAB
CUTTING EDGE TECHNOLOGY

VP Premium 2.0

The position of wear

After fitting in a frame, the lens is very often tilted so the vision becomes blurry. This is because it is not optimized accurately according to the position of wear.

Thanks to **OptiMesh** technology, **VP Premium 2.0** corrects the difference in power between the theoretical and the resulting prescription when the user is wearing a lens.

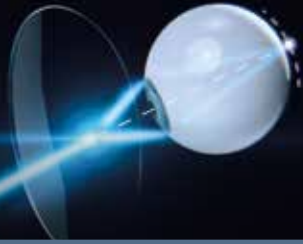

6 parameters are considered for the calculation



Fitting height Naso-pupillary distance Back Vertex distance



Wrap angle Pantoscopic angle Working distance

OptiMesh technology

When the eye rotates in a non-compensated lens surface, it appears oblique astigmatism.

To avoid that, OptiMesh technology uses Ray-Tracing technique to simulate how the lens behaves when it's placed in front of the eye of the wearer.

It also calculates what oblique astigmatism would be present in each point of the lens when the eye rotates.

VP Premium 2.0 offers excellence in visual acuity all over the lens surface.



With **VP Premium 2.0** you will notice visual sharpness in all gaze directions.

Personalized lens with differentiated features



Visual acuity in all gaze directions



More comfort while reading



Adapted to the frequent use of hand held devices



No restrictions on frame selection